




munevo DRIVE

**INNOVATIVE WHEELCHAIR CONTROL
VIA SMARTGLASS**

munevo DRIVE

Many people worldwide depend on an electric wheelchair. In addition, a large number of them are not able to control their wheelchair with their own hands.

munevo DRIVE is a Smartglass application that detects head movements with the help of Smartglass and translates them into steering signals. These signals are sent via Bluetooth wirelessly to the munevo adapter and from there to the control unit of the wheelchair. This enables hands-free control of the wheelchair by **munevo DRIVE**.

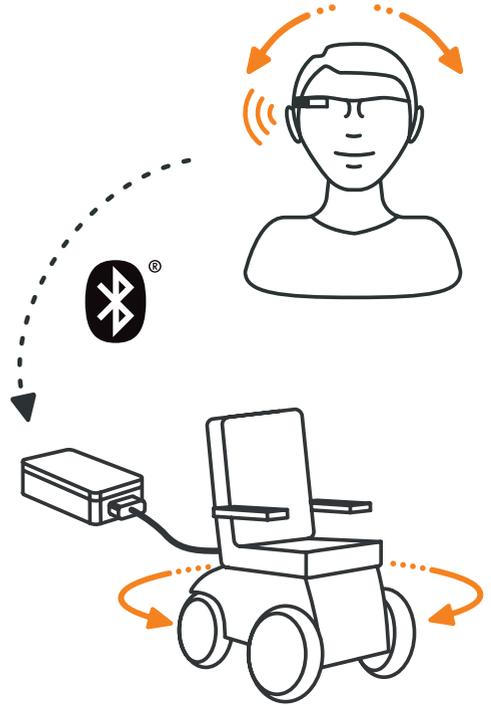

„An optimal solution so that I can operate the wheelchair myself despite my handicap.“



„The menu navigation is incredibly clear and very easy to learn.“

Functions



drive



take pictures



seat adjustment



share pictures



emergency message



Robot

Planned



Navigation



Smartphone



Smart-Home



Computer

Advantages - why **munevo DRIVE**?



Comfortable and modern

Due to its attractive design, the Smartglass controller is less stigmatizing than conventional solutions. The lightweight eyewear is comfortable to wear and very convenient for daily use.



Connectivity and more features

In the future, Bluetooth technology will make it possible to control a large number of other products wirelessly. For example, the user will be able to control his smart home or a robot arm via the **munevo App**; for a larger radius of action and more independence in everyday life.



Adaption

Unlike chin or occipital controls, the ability to control the wheelchair does not depend on the user's position in the wheelchair. The user can compensate for changes in the sitting position or the range of movement by calibrating the app.



Indoor and Outdoor

The Smartglass and the **munevo adapter** are waterproof according to IPX4 or IP65. Thus the system can be used both indoors and outdoors.



Safe

By integrating **Bluetooth** security, the system is protected against external access.



GDPR conformity

All user data is processed by **munevo** according to GDPR. Data is only recorded to improve the product or service.

Productdetails



Smartglass with frame:

- Google Glass Enterprise Edition
- Weight: 41 g (incl. frame)
- **munevo DRIVE** application pre-installed

munevo DRIVE adapter:

- Receiver of driving signals (via **Bluetooth**)
- Interface to wheelchair input module (Omni)



Contact

munevo GmbH

+49 (0)89 550649-04

info@munevo.com

www.munevo.com



munevo

Address:

munevo GmbH

Nordostpark 89

90411 Nürnberg

Postal address:

munevo GmbH

Agnes-Pockels-Bogen 1

80992 München



Brands:

munevo, **munevo DRIVE** and the **munevo** logos are registered trademarks of .

Bluetooth is a worldwide registered trademark of Bluetooth SIG, Inc.

Copyright © 2019 **munevo GmbH**